

Exp → Estimation of Iron by Potentiometric titration

### Aim

To estimate the amount of  $\text{Fe}^{2+}$  ion present in the given solution

### Apparatus required

Potentiometer assembly, 25 ml burette, 10 ml pipette, 250 ml beakers, standard flask, calomel and platinum electrodes.

### Reagents required

Ferrous ammonium sulphate,  
dil.  $\text{H}_2\text{SO}_4$ , std.  $\text{K}_2\text{Cr}_2\text{O}_7$

### Principle

Measurement of Emf of a electrochemical cell using potentiometer and change in Emf due to the chemical (redox) reaction is monitored.

In this potentiometric titration setup an indicator electrode (Pt electrode) & reference electrode (calomel) is coupled to form electrochemical cell for  $\text{Fe}^{2+}$  to  $\text{Fe}^{3+}$ .

$\text{Fe}^{2+}$  is oxidizing to  $\text{Fe}^{3+}$  as  $\text{K}_2\text{Cr}_2\text{O}_7$  is progressively added. Platinum electrode which is kept in contact with a mixture of  $\text{Fe}^{2+}$  &  $\text{Fe}^{3+}$  ions acts as a redox electrode (indicating redox rxn)


The reduction potential of this single electrode depends on the ratio of  $[Fe^{2+}/Fe^{3+}]$  initially. During the titration of  $Fe^{2+}$  in  $H_2SO_4$  medium with  $K_2Cr_2O_7$ , this ratio varies to a little extent at the beginning and suddenly near the end point. After the end point it the ratio changes very little. It can be noted that there is a sudden change in the ratio of  $[Fe^{3+}/Fe^{2+}]$  as the equivalence point is reached. This causes a sudden increase in the EMF of the cell at equivalence point.

### Procedure

- (I) The given ferrous iron solution is made upto a known volume (say 100 ml) in a SMF following the standard procedure with usual precautions.
- (II) Exactly 10 ml of the made up  $Fe^{2+}$  solution is pipetted out into a clean 100 ml beaker. About 10 ml of dilute  $H_2SO_4$  and 100 ml of distilled water are added to it.
- (III) A platinum electrode is dipped into this solution & coupled with the standard calomel electrode. The resultant cell is then incorporated into the potentiometric

Aim - To estimate the amount of  $\text{Fe}^{2+}$  ion present in the given solution.


The cell setup:

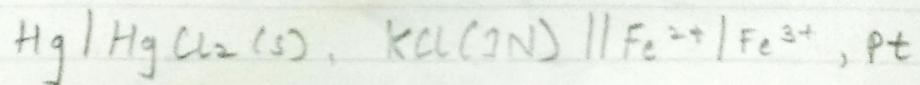

The chemical reaction:



Pilot titration

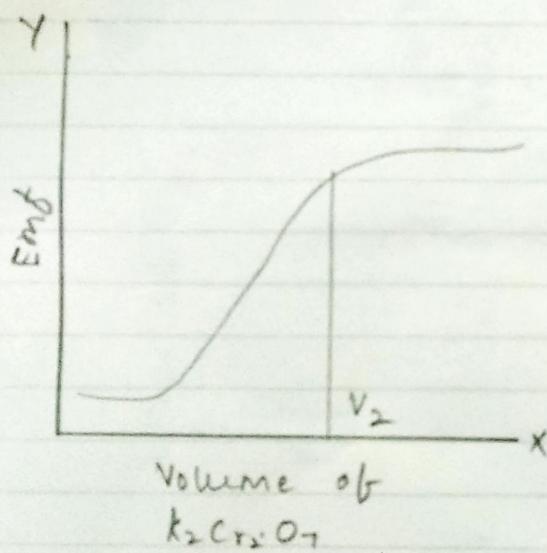

Fair titration

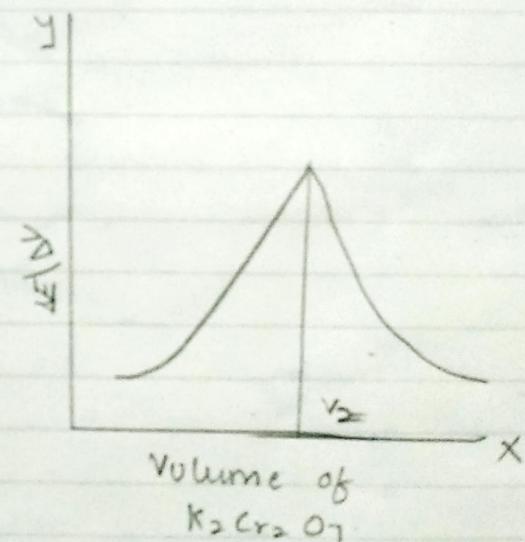


Table 1 : Pilot titration : FAS vs  $K_2Cr_2O_7$

FAS  $\rightarrow$  Ferrous Ammonium Sulphate

Sr. No.	Volume of $K_2Cr_2O_7$ (ml)	EMF (mV)	$\Delta E$ (mV)
1.	0	364	- Nit -
2.	1	380	16
3.	2	400	20
4.	3	423	23
5.	4	440	17
6.	5	460	20
7.	6	490	30
8.	7	530	40
9.	8	670	40
10.	9	690	20
11.	10.	700	11
12.	11	720	19
13.	12	740	20
14.	13	745	5
15.	14	760	15
16.	15	772	12

## Pilot titration

(Charvi Jain)  
RA2111047010113

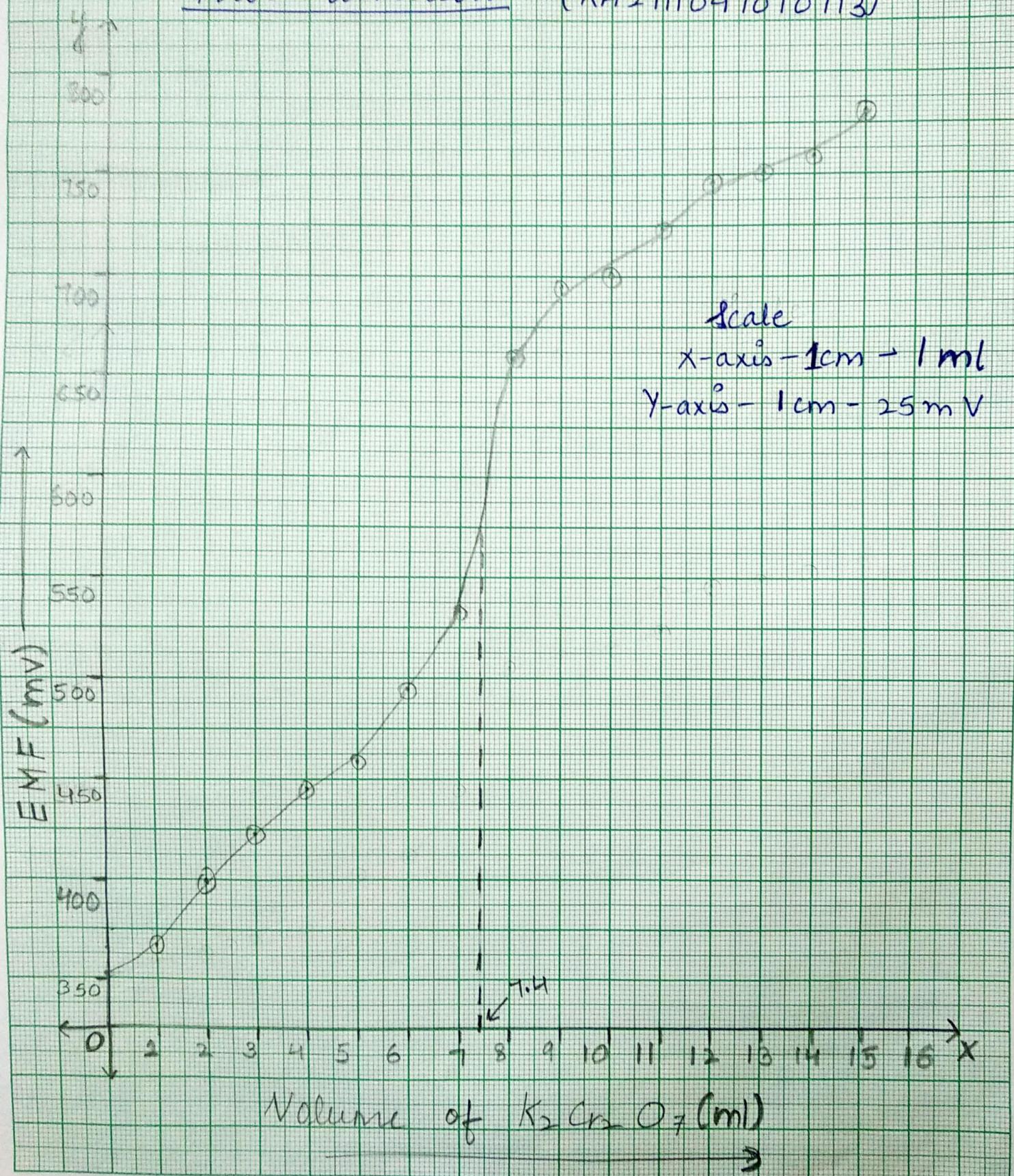


Table 2 : Fair Titration - FAS vs  $K_2Cr_2O_7$

Sr. No.	Volume of $K_2Cr_2O_7$ (ml)	EMF (mV)	$\Delta E$ (mV)	$\Delta V$ (ml)	$\Delta E / \Delta V$
1.	5	465	-	-	-
2.	5.2	470	5	0.2	25
3.	5.4	480	10	0.2	50
4.	5.6	491	11	0.2	55
5.	5.8	496	5	0.2	25
6.	6	500	4	0.2	20
7.	6.2	510	10	0.2	50
8.	6.4	521	11	0.2	55
9.	6.6	533	12	0.2	60
10.	6.8	543	10	0.2	50
11.	7	550	7	0.2	35
12.	7.2	560	10	0.2	50
13.	7.4	630	70	0.2	350
14.	7.6	650	20	0.2	100
15.	7.8	670	20	0.2	100
16.	8	680	10	0.2	50
17.	8.2	685	5	0.2	25
18.	8.4	690	5	0.2	25
19.	8.6	694	4	0.2	20
20.	8.8	704	10	0.2	50
21.	9	710	6	0.2	30

## Fair titration

(Charvi Jain  
RA2111047010113)

Y ↑

Scale

X-axis - 1 cm - 0.2 ml

Y-axis - 1 cm - 50 unit

300

250

200

150


100

50

$\Delta E / \Delta V$

0 5 5.2 5.4 5.6 5.8 6 6.2 6.4 6.6 6.8 7 7.2 7.4 7.6 7.8 8 8.2 8.4 8.6 8.8 9

Volume of  $K_2Cr_2O_7$  (ml) →

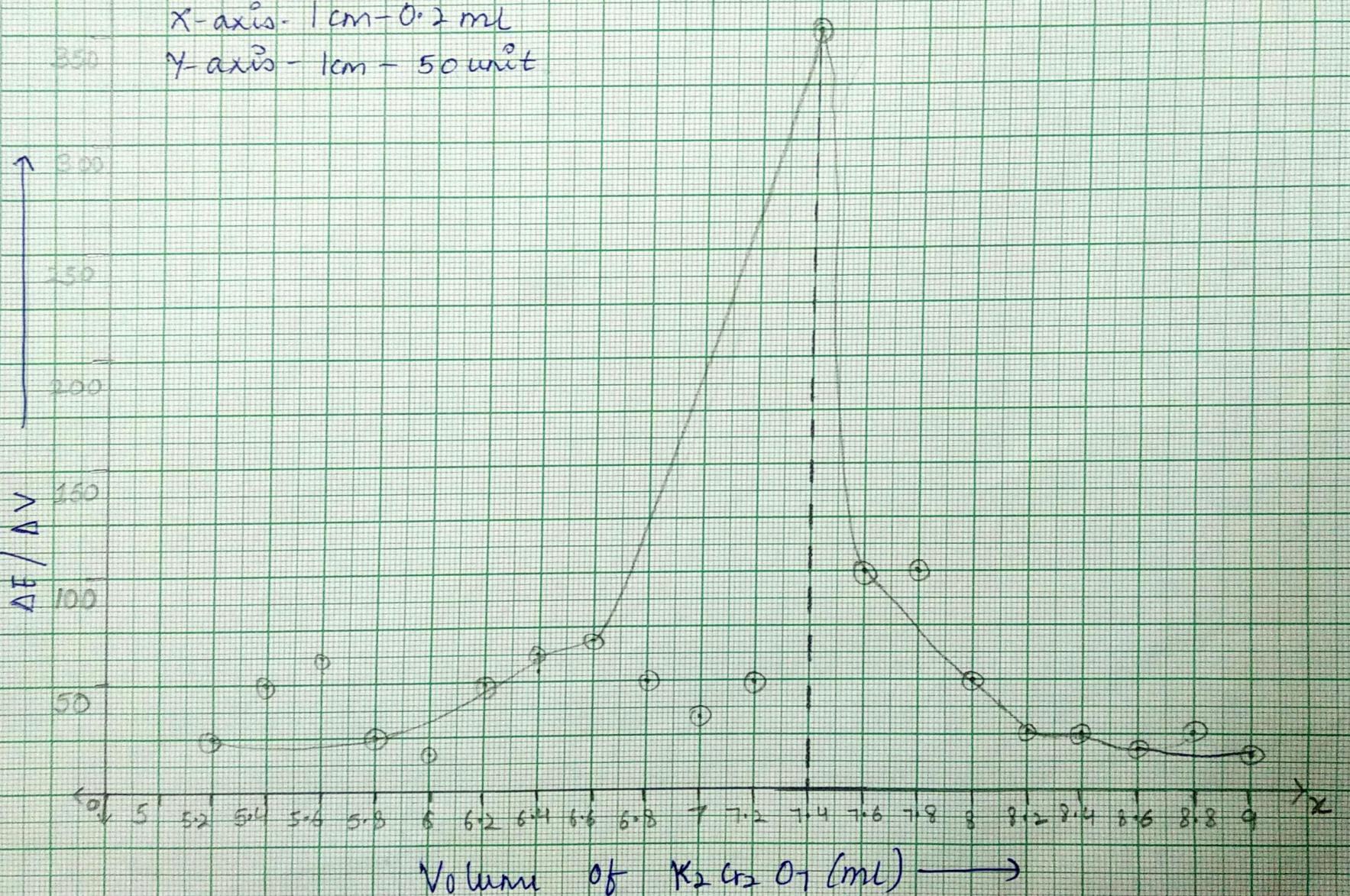

### \* Calculation:

$$V_1 = \text{Volume of pipette solution (FAS)} = 10 \text{ ml}$$

$$V_2 = \text{Volume of } K_2Cr_2O_7 = 7.41 \text{ ml}$$

$$N_2 = \text{Normality of } K_2Cr_2O_7 = 0.1 \text{ N}$$

$$\text{Normality of FAS (N)} = ?$$

$$\text{Strength of FAS} = N_1 = \frac{V_2 N_2}{V_1}$$
$$= \frac{7.41 \times 0.1}{10} = 0.074 \text{ N}$$

$$\boxed{\text{Normality} = 0.074 \text{ N}}$$

$$\text{Amount of } Fe^{2+} (\text{g/lit}) = \text{Eq.wt.} \times \text{Normality of } Fe^{2+}$$
$$= 55.85 \times 0.074 \text{ (g/lit)}$$

$$\boxed{\text{Amount of } Fe^{2+} = 4.132 \text{ g/lit}}$$

$$\text{Amount in } Fe^{2+} \text{ in 100 ml} = \text{Eq.wt.} \times \text{Normality of } Fe^{2+}$$
$$= 55.84 \times 0.074 = 0.0413 \text{ grams}$$

### \* Result


Weight of Iron in given solution = 0.0413 grams

circuit.

- (iv) Standard  $K_2Cr_2O_7$  solution, which is taken in a burette, is added in installments of 1ml into the beaker and the cell e.m.f. is measured after each addition by proper mixing.
- (v) The process is continued still or also well beyond the neutralization point as indicated by an abrupt change in the e.m.f.

The cell setup  $\rightarrow Hg | HgCl_2(s), KCl(N) || Fe^{2+} | Fe^{3+}, Pt$

The chemical reaction:

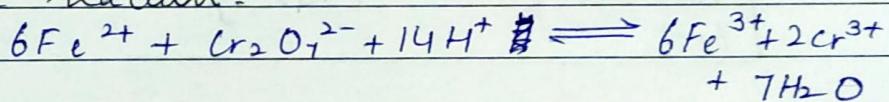

- \*(i) Note the volume  $K_2Cr_2O_7$  solution required for complete oxidation of  $Fe^{2+}$  solution from the plot of e.m.f. vs the volume of the standard  $K_2Cr_2O_7$  solution added. The range at which the endpoint lies may be evaluated.
- (ii) Calculate the normality of given  $Fe^{2+}$  solution using the formula  $N_1V_1 = N_2V_2$
- (iii) One more similar titration is performed, using 0.1 ml. portions of  $K_2Cr_2O_7$  (std.)

solution close to the end point (1ml on either side of the range) and tabulating the measured emf. corresponding to each addition.

- (IV) Graph is plotted (i) Emf vs Volume of  $K_2Cr_2O_7$  addition  
 (ii)  $\Delta E / \Delta V$  vs volume of  $K_2Cr_2O_7$  addition
- (V) The exact end point can be determined from the plot of  $\Delta E / \Delta V$  versus volume of standard  $K_2Cr_2O_7$  solution.